
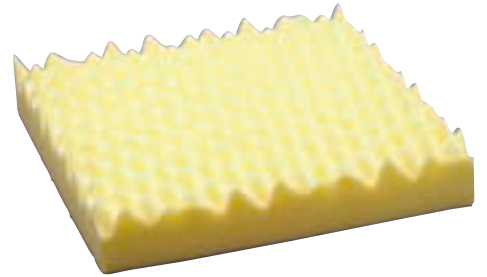


Cushions

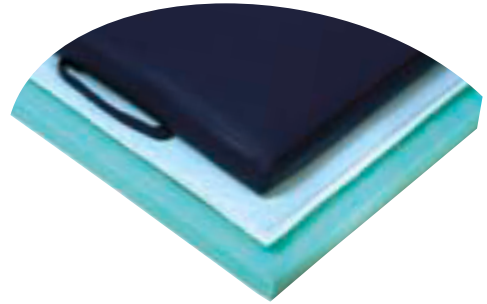

GELIDIS CUSHIONS

In compliance with new antidecubitus guidelines, the *GELIDIS* cushion range offers all the required features and optimal comfort of use.



Our pillows are composed of a viscoelastic gel module and a base of three pre-formed layers of high resilience foam, one of which is anatomically shaped at the ischion level. The adequate gel thickness preserves its stiffness and its excellent elasticity ensures optimal load distribution, thus actively contributing to bedsores prevention. They are covered with a micropore breathing cover and non-slip fabric on one side, easily removable for maintenance.



Washing recommendations

- High resilience foam and gel: soapy water
- Jersane cover : 60°C
 - Somy cover: 90°C
 - Fur cover : 60° C

- 03/2553H JERSANE cover for GELIDIS PLUS
- 03/2554H SOMY cover for GELIDIS
- 03/2557H SOMY cover for pierced GELIDIS
- 03/2558H Fur cover for GELIDIS

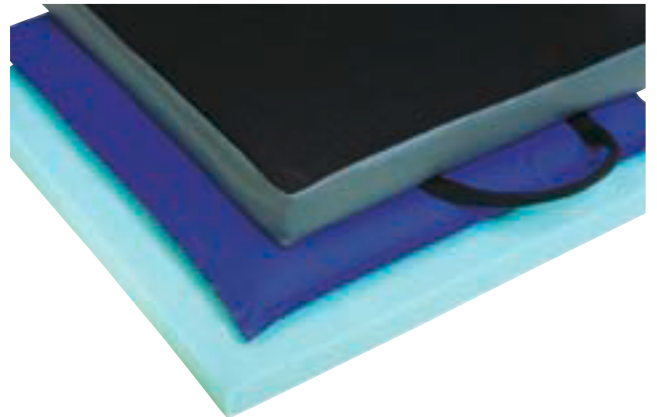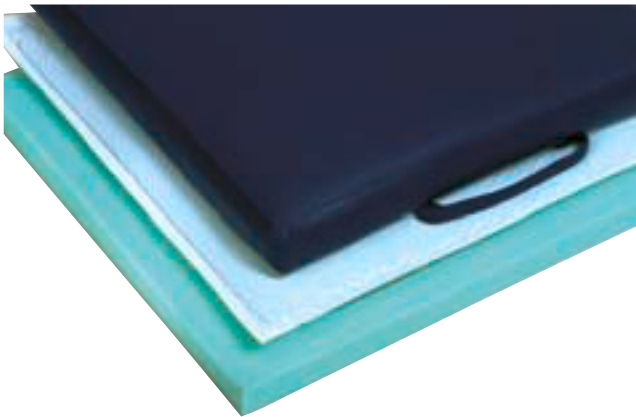

GELIDIS and pierced GELIDIS : supplied with SOMY breathing cover, polyurethane-lined, bacterium-resistant fabric (Bio-Pruf treatment), M1 fire retardant classified. Washable at 90°C.

- GELIDIS : 40 x 40 x 5 cm • Weight : 3,5 kg
- Pierced GELIDIS : 40 x 40 x 3,5 cm • Weight : 3,2 kg

GELIDIS PLUS : Supplied with a JERSANE cover, combining cotton resistance with elastane elasticity.

- M1 fire retardant classified. Washable at 60°C.
- GELIDIS PLUS : 40 x 40 x 6 cm • Weight : 4,7 kg



03/2553
GELIDIS PLUS



03/2554
GELIDIS



03/2557
PIERCED GELIDIS